



**HUMAN CARE**  
Your Life. Your Way.



**Mobility Aids**

# Introduction

We produce and provide mobility solutions for people with special needs. Our products are made to serve as a natural part of life to all our users.

Human Care has a proud history in patient care. Founded in 1993, our twenty plus years of experience is expressed through a broad selection of premium products offering high customer value and satisfying all customer needs.

Our main product areas are: Mobility Aids, Lifting Solutions, Healthcare Beds and Convertible Chairs for people with reduced mobility.

With headquarters in Stockholm, Sweden, and subsidiaries in the USA, Canada, the Netherlands and Australia, we collaborate with distributors to provide solutions in more than 40 countries.

Our promise is safety, quality and usability. In all of our products, we exercise only the highest quality to make sure that every detail, every aspect, every function will deliver that promise made to our users.

**Welcome to Human Care.**

## Content

### Rollators

Leia	8
neXus 3	12
neXus 1, 2	15
Carl-Oskar	16
Rebel	18
Parade	21
Fellow Supreme	22
Fellow Classic	22

---

### RA Rollators and Trolley Walkers

Carl-Oskar RA	26
Rebel RA	26
Staffan RA	26
Carl-Oskar TW	28
Rebel TW	28
Staffan TW	28
Svea TW	28

## Your life. Your way.

That is what we have been offering people with reduced mobility for over 20 years. Your Life. Your Way. Our focus is simple. We work for our users. Our products serve to be a natural part of life to anyone who might need them. Meeting that need is what drives and motivates us. Our mission is to improve the freedom of movement for people with special needs. Our solutions provide reliability, high quality, innovation and state-of-the-art technology. In every detail, in every aspect, our products deliver what they promise.

“Our mission is to improve the quality of care and freedom of movement for people with special needs”



## Benefits with rollators

The rollator is a mobility aid that greatly enhances the user's autonomy. The rollator allows people of reduced strength and balance a social life – as well as exercise, improving blood circulation and general health. The rollator gives support and relief during a walk – as well as serving as a seat when the user needs a rest.


By allowing the users to decide themselves when to go out or to socialize, the rollator greatly enhances the users' self-esteem and quality of life. Thus, the rollator – directly and specifically – contributes to the strengthening of public health by allowing people to live life their way.

Today, our product range encompasses more than 30 different rollator models: indoor, outdoor, trolleys, RA walking tables – as well as a wide range of accessories that meet all user demands.

## Care and Quality

A user of a product from Human Care is a person who feels very confident with the product's quality and reliability. At all times. Every step of the way. The reason behind this is our profound commitment to quality, safety and usability.

Firstly, our products are designed and manufactured to meet the requirement of relevant ISO and EN standards. Secondly, every single product is given extensive attention at every stage of the design process. Thirdly, in developing our products, our design team collaborates with globally leading physiotherapists, creating state-of-the-art products which all share the same ultimate goal: the benefit and enjoyment of their users. Allowing them to live life - their way, every step of the way.

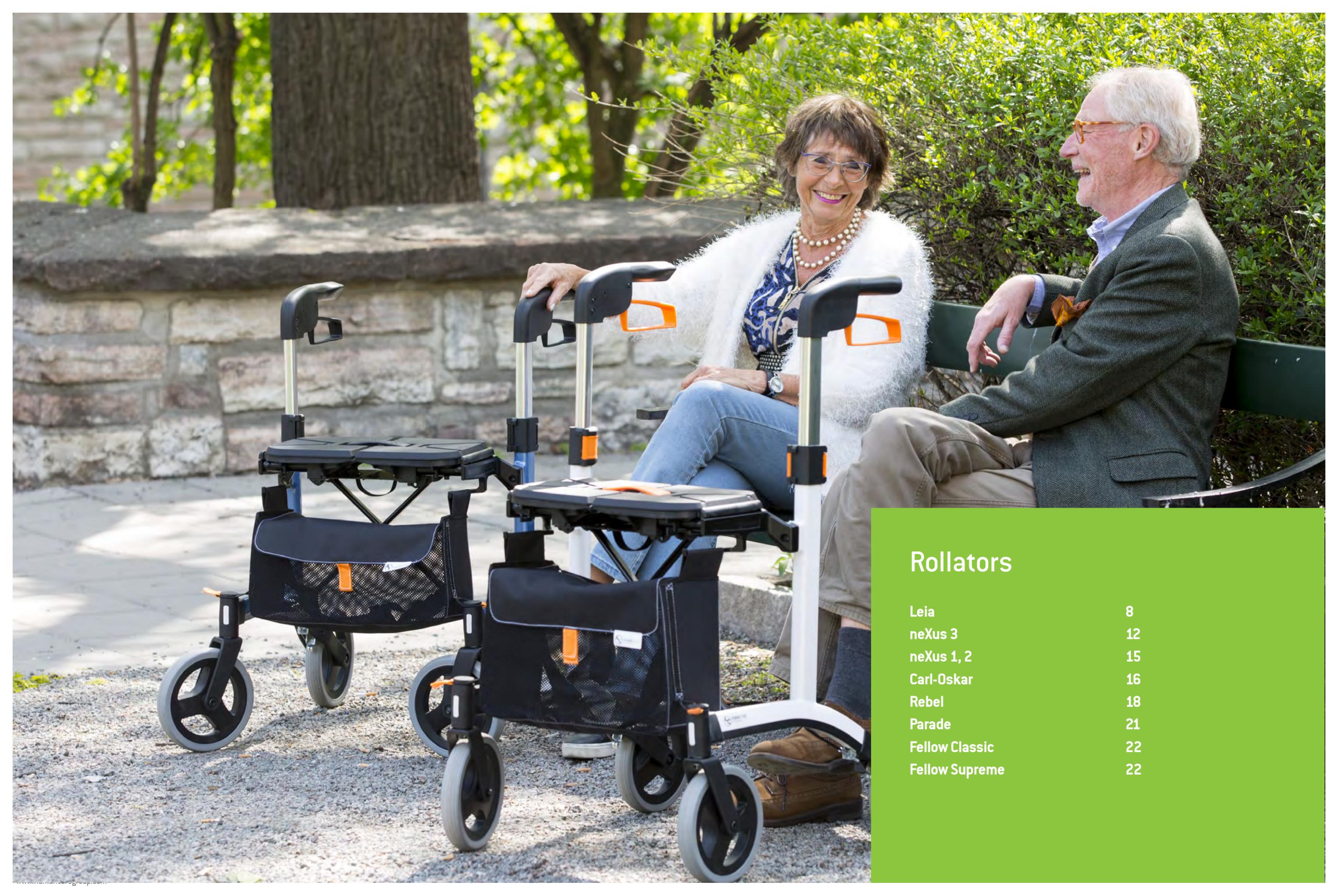

## Rollators

Leia	8
neXus 3	12
neXus 1, 2	15
Carl-Oskar	16
Rebel	18
Parade	21
Fellow Classic	22
Fellow Supreme	22

# Leia

The Leia is an exciting new addition to our family of rollators. A modern series of rollators, the Leia is ideal for those who place high demands on both design, quality and function. These rollators have been designed with a focus on design and usability and are available in three different colors and three seat heights to allow you to find something ideal for you.

Thanks to the patented cross folding that enables ample walking space, the new Leia offers the best possible ambulatory pattern ensuring less stress on both back as well as shoulders. The Leia is lightweight and is very easy to maneuver yet has maximum stability.



Leia white

Leia grey

Leia blue

## Superior walking ergonomics and maneuverability

With Leia's unique frame design, the seat is placed further forward into the front of the rollator. This enables an ample walking space and the best possible ambulatory pattern which is easy on both the back as well as the shoulders. The Leia is lightweight and very easy to maneuver yet still offers maximum stability. The sturdy yet flexible frame design helps to absorb bounces or fluctuations from ground impact.

## Excellent design and comfortable resting spot

The Leia has been designed with focus on design, quality and usability and is available in three different colours and three heights to allow you to find something ideal for you. The contoured, cushioned seat offers a very comfortable resting place. For extra comfort and support while resting, a padded back rest is also available as an accessory<sup>1</sup>.

## Lightweight and easy to fold and transport

The patented cross-folding design makes the rollator very easy to fold and store and thanks to the low weight it is also easy to lift and transport. It also stands independently when folded and takes up minimal storage space. The soft bag comes as a standard option and folds together with the rollator.

## Accessories

	<b>Tray</b> Practical for resting light weight items. Max load: 2 kg. Part no: 89747		<b>Cane holder</b> Simplifies the transportation of a cane or umbrella. Part no: 80603		<b>Backrest</b> Gives support for the back while sitting. Wide design for extra support. Part no: 89746
	<b>Slow down brake</b> Prevents the rollator from moving too fast and enables a controlled walking speed for the user. Part no: 80620		<b>Accessory bag</b> Attaches to the back side of the seat pointing towards the user. Ideal for storing personal belongings. Part no: 89654		<b>Curb climber</b> Helps raise the rollator onto a curb. Part no: left: 80612, right: 80613

## Technical Specifications

Description	Part no. global	Part no Canada	Handle Height Min.	Handle Height Max.	Seat height	Width	Width between handles	Outside Handle Width	Wheel size	Weight	Max user weight
Leia high white	70949	90259	80 cm / 31.5"	97 cm / 38.2"	62 cm / 24.4"	61 cm / 24.1"	44 cm / 17.3"	50 cm / 19.6"	20x3 cm / 7.9x1.2"	7,9 kg / 17.4 lbs	150 kg / 330 lbs
Leia high blue	89729	90260									
Leia high grey	89720	90261									
Leia medium white	89762	90267	71 cm / 28"	88 cm / 34.7"	57 cm / 22.4"	61 cm / 24.1"	44 cm / 17.3"	50 cm / 19.6"	20x3 cm / 7.9x1.2"	7,6 kg / 16.8 lbs	150 kg / 330 lbs
Leia medium blue	89814	90266									
Leia medium grey	89821	90265									
Leia low white	89795	90262	66 cm / 26"	76 cm / 30.0"	52 cm / 20.5"	61 cm / 24.1"	44 cm / 17.3"	50 cm / 19.6"	20x3 cm / 7.9x1.2"	7,5 kg / 16.5 lbs	150 kg / 330 lbs
Leia low blue	89833	90264									
Leia low grey	89840	90263									

Note: Seat sizing reference: High is the equivalent of Standard, Medium the equivalent of Low, Low the equivalent of Super Low. <sup>1</sup>Backrest is included in North America.



### The importance of proper walking

Walking involves several body parts all interacting together - your head carriage, shoulder girdle, arm swing, hip movement, knee action and foot placement all work together while walking. Any disconnect in this interaction, whether a result of disability or age, will require some type of assistance to ensure proper walking is maintained.

Studies have shown that loss of function as a result of disease or age impairs sensorimotor functions and changes both gait velocity and pattern compromising overall mobility. As a result, a compensatory walking pattern can occur causing undue stress on the back and shoulders as well as additional encumbrances on other muscle groups.

### Using a rollator that encourages proper walking patterns

Being able to “step into” a rollator ensures the rollator is not too far away from the user while taking steps. Traditional rollators the user tends to bend over to reach the handles. With a seat that does not obstruct the gait pattern and accommodates the wider and shorter step lengths taken, it allows for a faster and more normal gait pattern, mimicking the correct pattern that should be used.



# neXus 3

The best-selling neXus 3 is an innovative, high quality rollator. The neXus 3 features an "X" frame design that folds from side to side which makes it quick and easy to fold. It also allows the rollator to stand independently when folded and ensures that it requires minimal space.

The cable-free brakes ensure that no cables can be caught on objects around the home, thus increasing safety. The rollator has a padded back strap that can flip to allow for two-position seating and a double bearing front fork assembly for smooth rolling. The neXus 3 comes with a soft bag that can remain on the rollator even when folded. Optional kits that include a colored soft bag with a coordinating back strap and handle are also available to personalize your rollator!



## Cable-free braking system

Thanks to the cable-free braking system there is less risk of injury and at the same time increased safety since no cables catch on objects around the home.

## Superior walking ergonomics

The cross-folding design enables an ample walking space which results in the best possible walking pattern for the user. The padded, contoured, cushioned seat offers a very comfortable resting spot. For extra comfort and support a padded back rest comes as standard feature on the rollator.

## Lightweight and easy to fold even with the bag on the rollator

The cross-folding design makes the rollator very easy to fold and store. It also stands independently when folded. The soft bag comes as a standard option and folds together while remaining on the rollator. It has convenient pockets for storing items such as a tablet or mobile phone.



Accessory kit purple

Accessory kit tartan

Accessory kit blue

## Accessories

	<b>Tray</b> Practical for resting light weight items. Part no: 80702		<b>Hard basket</b> Basket ideal for shopping. Max capacity: 10 lbs. Part no: 4510		<b>Soft bag</b> Soft bag that folds with the rollator. Included with product or ordered. Max capacity: 10 lbs. Part no. 8050105 (5 pieces).
	<b>Slow down brake</b> Prevents the rollator from moving too fast and enables a controlled walking speed for the user. Part no: 4950		<b>Backrest</b> Gives support for the back when sitting. Included with product.		<b>Accessory kits</b> Accessory kit of soft bag, back strap and handle. Purple, part no. 89991 Tartan, part no. 89992 Blue, part no: 89990
	<b>Cane holder</b> Simplifies the transportation of a cane or umbrella. Part no: 80703		<b>Curb climber</b> Allows rollator to be slightly lifted onto a curb. Included with product.		

## Technical Specifications

Description	Part no.	Handle Height min	Handle Height max	Seat height <sup>1</sup>	Width	Width between handles	Wheel size	Weight	Max user weight
neXus 3 super low with soft bag	460018	Reg: 70 cm/27.5" Tall: 78 cm/30.5"	Reg: 78cm/30.5" Tall: 85cm/33.5"	47 cm / 18.5"	61 cm / 24"	45 cm/17.5"	20 x3 cm /7.9"x1.2"	6.8 kg / 15 lbs	150 kg / 330 lb
neXus 3 Low with soft bag	460022	Reg: 75 cm/29.5" Tall: 78 cm/ 30.5"	Reg: 85 cm/33.5" Tall: 93 cm/36.5"	53 cm / 21"	61 cm / 24"	45 cm/17.5"	20 x3 cm /7.9"x1.2"	6.8 kg / 15 lbs	150 kg / 330 lb
neXus 3 Standard with soft bag	460024	Reg: 83 cm/ 32.5" Tall: 83 cm/ 32.5"	Reg: 93cm/36.5" Tall: 100 cm/39.5"	61cm / 24"	61 cm / 24"	45 cm/17.5"	20 x3 cm /7.9"x1.2"	7.3 kg / 16 lbs	150 kg / 330 lb

<sup>1</sup>Add 1" for cushion, For tall handles, add "T" at end of product code.



# neXus 1, neXus 2

## neXus 1

The nexus 1 walking aid features 6"/15cm wheels that are ideal for indoor use. Its ultra-compact "X" folding system allows it to stand independently when folded. The unique design allows the user to walk with the rollator partially folded, making it ideal for those smaller indoor spaces.



neXus 1

## neXus 2

The neXus 2 is a lightweight rollator with great stability. It features an "X" frame design that folds from side to side which makes it quick and easy to fold. It also allows the rollator to stand independently when folded and ensures that it requires minimal space.



neXus 2

## Accessories

	<b>Tray</b> Practical for resting light weight items. Part no: 80702		<b>Hard basket</b> Basket ideal for shopping. Part no: 4510. Included with rollator.		<b>Curb climber</b> Helps lift the rollator onto a curb. For neXus 2 only. Part no. 4878
	<b>Cane holder</b> Attaches quickly to hold standard canes. Part no: 80703		<b>Backrest</b> Gives support for the back while sitting. Included with product. Part no: P49050		<b>Soft bag</b> Soft bag with convenient pockets, folds with the rollator. Part no. 80501Q5 (5 pieces).

Note: Nexus accessory kits can be used on the neXus 1 or neXus 2 with purchase of clamp. Part no: 90276.

Technical Specifications									
Description	Part no.	Handle height min	Handle height max	Seat height <sup>1</sup>	Width	Width between handles	Wheel size	Weight	Max user weight
neXus 1, Super Low	452018	71 cm / 28"	81 cm / 32"	45 cm / 17.5"	64 cm / 25"	45 cm / 17.5"	15 cm x 3 cm / 6" x 1.2"	5.9 kg / 13 lbs	113.4 kg / 250 lbs
neXus 1, Low	452022	78 cm / 30.5"	88 cm / 34.5"	51 cm / 20"				5.9 kg / 13 lbs	
neXus 1, Standard	452024	85 cm / 33.5"	95 cm / 37.5"	58 cm / 23"				6.3 kg / 14 lbs	
neXus 2, Super Low	451518	74 cm / 29"	84 cm / 33"	47 cm / 18.5"	64 cm / 25"	45 cm / 17.5"	20x3 cm / 7.9"x1.2"	6.8 kg / 15 lbs	113.4 kg / 250 lbs
neXus 2, Low	451522	80 cm / 31.5"	90 cm / 35.5"	53 cm / 21"				6.8 kg / 15 lbs	
neXus 2, Standard	451524	88 cm / 34.5"	98 cm / 38.5"	61 cm / 24"				7.3kg / 16 lbs	

Only distributed within Canada. <sup>1</sup> add 1" for seat pad]

# Carl-Oskar

The Carl-Oskar is one of Human Care's best known and top-selling rollators with excellent stability. It is a safe and reliable all-round rollator available in four different widths and with a weight capacity of 150 kg. Its key feature is the flexibility which allows it to be stable even on rough or uneven surfaces. With large 11" wheels, the Carl-Oskar is ideal for rural areas. The Carl-Oskar is also easy to maneuver and has ergonomic handles for a stable and comfortable grip.

Folding is easy and the rollator stands independently in the folded position. The Carl-Oskar is also available as a reinforced model with a maximum weight capacity of 275 kg/606.3 lbs.



## Outstanding stability


The large stable wheels make the rollator roll very easily yet remain stable on virtually any surface. The flexibility in the chassis allows the rollator to be stable even on rough or uneven surfaces. This makes the rollator one of the most versatile rollators on the market.

## Easy to fold and store

The Carl-Oskar rollator folds very easy from side to side and can stand securely and independently when folded. The rollator also takes little space when folded and is very practical to transport.


## Ergonomic handles and easy braking

In order to offer you the best ergonomics, the handles feature an ergonomic design that ensures even distribution of weight for a stable, comfortable grip. In addition the handles are angle-adjustable, so that you can find the optimal position for your hands.

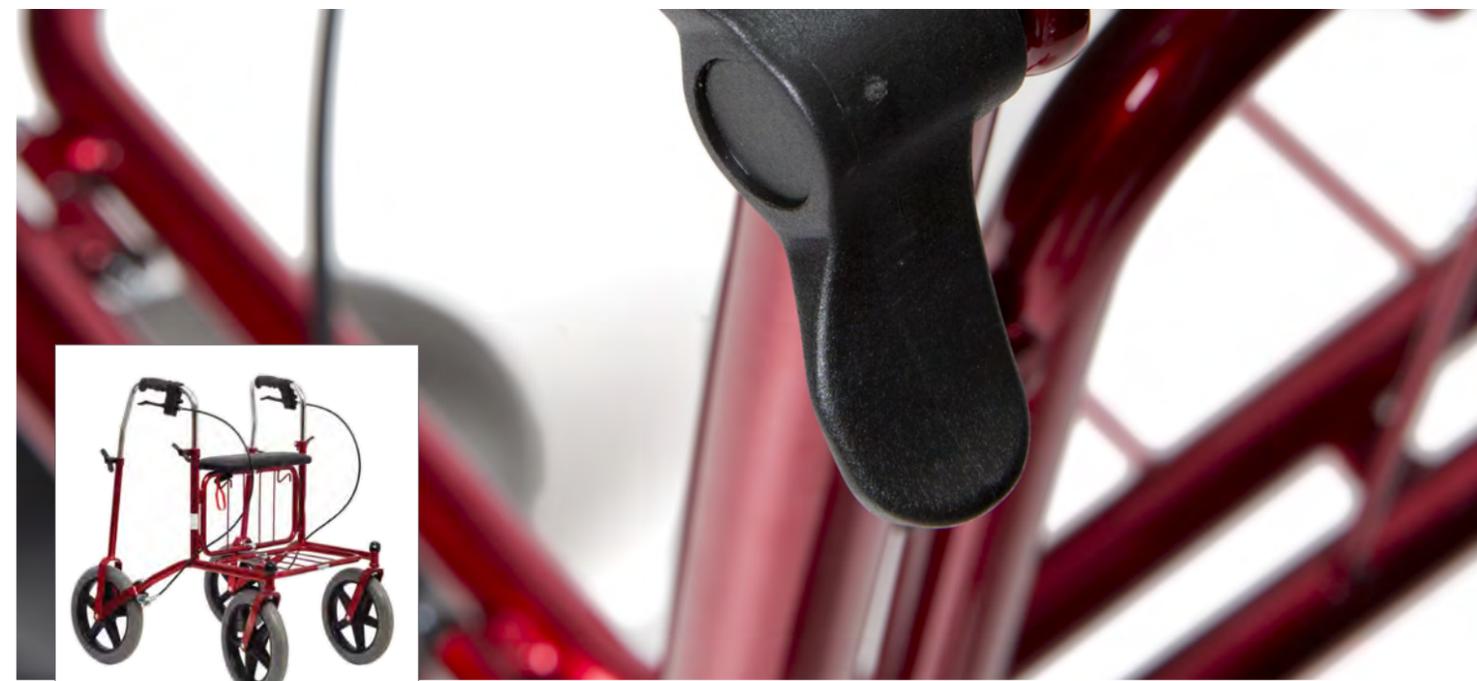

Carl-Oskar reinforced

## Accessories

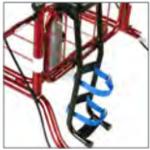	<b>Tray</b> Practical for resting light weight items. Part no. for 57 cm models: 10651, Part no. for 62cm, 67cm and 72cm models: 10602		<b>One-hand brake</b> Allows the user to control brake functions with one hand. Part no: 10652		<b>Splash guard</b> Protects legs and clothes. Sold in pairs. Part no: 10608
	<b>Slow down brake</b> Prevents the rollator from moving too fast and enables a controlled walking speed for the user. Part no: 10620		<b>Backrest</b> Supports the back while resting. Part no. for 57 cm, 62cm and 67 cm models: 10621. Part no. for 72cm model: 10655		<b>Oxygen holder</b> Holds oxygen tube firmly in place. Part no: 10615
	<b>Cane holder</b> Simplifies the transportation of a cane. Part no: 10603		<b>Curb climber</b> Helps raise the rollator onto a curb. Part no: 10612		<b>Basket with handle</b> Basket ideal for shopping. Part no: 10601


## Technical Specifications

Description	Part no.	Handle height min	Handle height max	Seat height	Width	Width between handles	Wheel size	Weight	Max user weight
Carl-Oskar 57cm/22.4"	11150	77cm	95 cm	64 cm	57 cm / 22.4"	37 cm / 14.6"	30 x 5 cm / 11.8"x1.8"	11.1 kg / 24.5 lbs	150 kg / 330.7 lbs
Carl-Oskar 62cm/24.4"	11250	77cm / 30.3"	95 cm / 37.4"	64 cm / 25.2"	62 cm / 24.4"	42 cm / 16.5"	30 x 5 cm / 11.8"x1.8"	11.3 kg / 24.9 lbs	330.7 lbs
Carl-Oskar 67cm/26.4"	11350				67 cm / 26.4"	47 cm / 18.5"		11.4 kg / 25.1 lbs	
Carl-Oskar 72cm/28.3"	11450				72 cm / 28.3"	52 cm / 20.5"		11.5 kg / 25.3 lbs	
Carl-Oskar inflatable tires, 57cm/22.4"	11151	77cm	95 cm	64 cm	57 cm / 22.4"	37 cm / 14.6"	30 x 5 cm / 11.8"x1.8"	11.1 kg / 24.5 lbs	150 kg / 330.7 lbs
Carl-Oskar inflatable tires, 62cm/24.4"	11251	77cm / 30.3"	95 cm / 37.4"	64 cm / 25.2"	62 cm / 24.4"	42 cm / 16.5"	30 x 5 cm / 11.8"x1.8"	11.3 kg / 24.9 lbs	330.7 lbs
Carl-Oskar inflatable tires, 67cm/26.4"	11351				67 cm / 26.4"	47 cm / 18.5"		11.4 kg / 25.1 lbs	
Carl-Oskar inflatable tires, 72cm/28.3"	11451				72 cm / 28.3"	52 cm / 20.5"		11.5 kg / 25.3 lbs	
Carl-Oskar reinforced, 62cm/24.4"	61225	77cm	95 cm	64 cm	62 cm / 24.4"	42 cm / 16.5"	30 x 5 cm / 11.8"x1.8"	13.2 kg / 29.1 lbs	275 kg / 606.3 lbs
Carl-Oskar reinforced, 67cm/26.4"	61325	77cm / 30.3"	95 cm / 37.4"	64 cm / 25.2"	67 cm / 26.4"	47 cm / 18.5"	30 x 5 cm / 11.8"x1.8"	13.3 kg / 29.3 lbs	606.3 lbs
Carl-Oskar reinforced, 72cm/28.3"	61425				72 cm / 28.3"	52 cm / 20.5"		13.4 kg / 29.5 lbs	

# Rebel

The high quality Rebel rollator is compact and easy to use both indoors or outdoors. The Rebel is safe and reliable and available in four different widths to accommodate a variety of users. With a weight capacity of 140kg, it is stable and easy to maneuver even on rough or uneven surfaces. The Rebel features comfortable, ergonomic handles for a stable and comfortable grip.

Folding is easy and the rollator stands independently in the folded position. The Rebel is also available in a reinforced model with a maximum user weight capacity of 250kg.

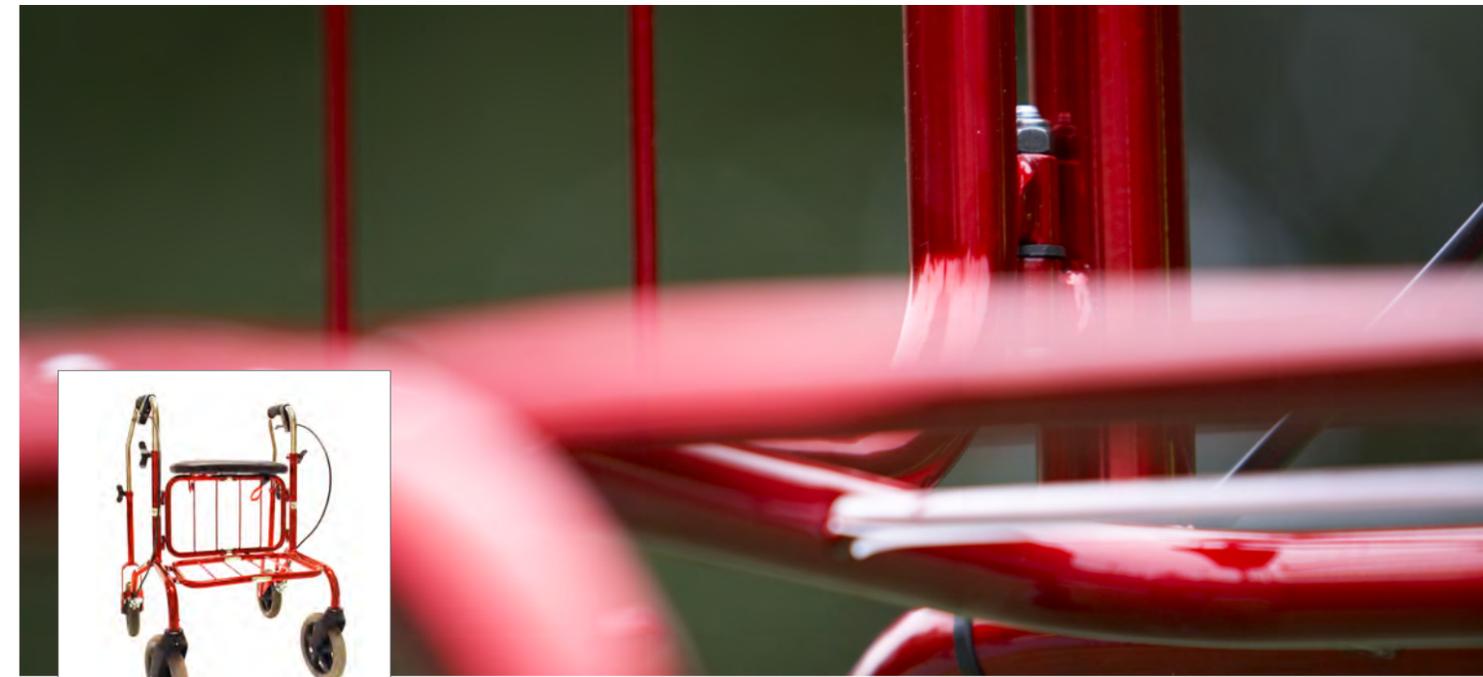

Rebel reinforced

## Perfect both indoors and outdoors

At home or on the go, the Rebel rollator is the perfect choice! Its flexible frame combined with stable wheels allow the rollator to move smoothly on virtually any surface.


## Compact design

The Rebel features a side folding system that can stand securely and independently while folded, making it very easy to transport.

## Ergonomic handles and easy braking

In order to offer you the best ergonomics, the handles offer a stable and comfortable grip. They are ergonomic and distribute the weight you put on them evenly. In addition the handles are angle-adjustable, so that you can find your optimal hand position.

## Accessories

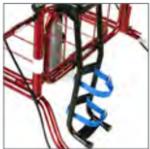	<b>Tray</b> Practical for resting light weight items. Part no. for 57 cm models: 10651, Part no. for 62cm, 67cm and 72cm models: 10602		<b>One-hand brake</b> Allows the user to control brake functions with one hand. Part no: 10652		<b>Oxygen holder</b> Holds oxygen tube firmly in place. Part no: 10615
	<b>Slow down brake</b> Prevents the rollator from moving too fast and enables a controlled walking speed for the user. Part no: 10620		<b>Backrest</b> Supports the back when resting. Part no. for 57 cm, 62cm and 67 cm models: 10621. Part no. for 72cm model: 10655		<b>Basket with handle</b> Basket ideal for shopping. Part no: 10601
	<b>Cane holder</b> Simplifies the transportation of a cane. Part no: 10603		<b>Curb climber</b> Helps raise the rollator onto a curb. Part no: 10612		

Technical Specifications									
Description	Part no.	Handle height min	Handle height max	Seat height	Width	Width between handles	Wheel size	Weight	Max user weight
Rebel 57cm/22.4"	22140	75 cm	930 mm	62 cm / 24.4"	57 cm / 22.4"	37 cm / 14.6"	20 x 4 cm / 7.9" x 1.6"	9.5 kg / 20.9 lbs	140 kg / 309 lbs
Rebel 62cm/24.4"	22240	/29.5"	/36.6"		62 cm / 24.4"	42 cm / 16.5"		9.6 kg / 21.2 lbs	
Rebel 67cm/26.4"	22340				67 cm / 26.4"	47 cm / 18.5"		9.7 kg / 21.4 lbs	
Rebel 72cm/28.3"	22440				72 cm / 28.3"	52 cm / 20.5"		9.9 kg / 21.8 lbs	
Rebel low, 57cm/22.4"	22141	68 cm	860 mm	54 cm / 21.3"	57 cm / 22.4"	37 cm / 14.6"	20 x 4 cm / 7.9" x 1.6"	9.2 kg / 20.3 lbs	140 kg / 309 lbs
Rebel low, 62cm/24.4"	22241	/26.8"	/33.8"		62 cm / 24.4"	42 cm / 16.5"		9.4 kg / 20.7 lbs	
Rebel low, 67cm/26.4"	22341				67 cm / 26.4"	47 cm / 18.5"		9.5 kg / 20.9 lbs	
Rebel low, 72cm/28.3"	22441				72 cm / 28.3"	52 cm / 20.5"		9.7 kg / 21.4 lbs	
Rebel reinforced, 62cm/24.4"	71221	75 cm	930 mm	62 cm / 24.4"	62 cm / 24.4"	47 cm / 18.5"	20 x 4 cm / 7.9" x 1.6"	11.5 kg / 25.4 lbs	250 kg / 551 lbs
Rebel reinforced, 67cm/26.4"	71321	/29.5"	/36.6"		67 cm / 26.4"	52 cm / 20.5"		11.6 kg / 25.6 lbs	
Rebel reinforced, 72cm/28.4"	71421				72 cm / 28.3"	52 cm / 20.5"		11.8 kg / 26.0 lbs	

## Parade

The Parade is designed for easy manoeuvrability even in small spaces. Despite its lightweight frame, it is still very stable and offers outstanding value for money. The maximum user weight capacity is 140 kg.

The Parade has ergonomic handles that ensure a comfortable and stable grip. It is easily folded and the rollator stands without support in the folded position.



### Designed for easy manoeuvrability

The Parade rollator is the perfect choice when looking for a lightweight rollator. Its design ensures easy manoeuvrability, even in the tightest spaces. Despite the lightweight frame, the Parade is still stable even on rough or uneven surfaces. The ergonomic handles ensure a comfortable and stable grip and the walker easily folds and stands independently in this position.

### Accessories



**Tray**  
Practical for resting light weight items.  
Part no. 10664



**Backrest**  
Supports the back when resting.  
Part no. 10621



**Basket**  
Basket ideal for shopping.  
Part no: 10662



**Cane holder**  
Simplifies the transportation of a cane.  
Part no: 10675



**Curb climber**  
Helps raise the rollator onto a curb.  
Part no: 10661

### Technical Specifications

Description	Part no.	Handle Height min	Handle Height max	Seat height	Width	Width between handles	Weight	Wheel size	Weight capacity
Parade	85500	74 cm / 29.1"	95 cm / 37.4"	60 cm / 23.4"	59 cm / 23.1"	43 cm / 16.9"	8.8 kg / 19.4 lbs	20 x 4 cm (7.8"x1.6")	140 kg / 309 lbs

# Fellow Series

## Fellow Supreme

The Fellow Supreme combines high quality ergonomic design with a low centre of gravity for increased safety and stability. Its vast handle height adjustment range and high weight capacity of 150kg make it easily suitable for most users.

With its unique fork construction, the Fellow Supreme is easy to maneuver. The soft braking system ensures that minimal hand pressure is required, allowing users with reduced hand functions to easily lock the wheels in place. The Fellow Supreme is also very easy to fold with a one-hand operation and can stand independently when folded.



Fellow Supreme

## Fellow Classic

The Fellow Classic is a high quality rollator ideal for both indoor and outdoor use. Safe and reliable it is available in five models: standard, low, high, wide and narrow.

The Fellow Classic's unique fork construction allows for easy manoeuvrability. Like the Supreme, it features a soft braking system ideal for those users with limited hand mobility and an easy-to-use folding system that will accommodate a variety of users. Weight capacity 150kg.



Fellow Classic

## Superior walking ergonomics

The forward facing handles allow the user to walk more upright which is easy on both back as well as shoulders. The Fellow series of rollators are also very easy to maneuver while also offering maximum stability. The sturdy seat flips up for added walking space or access to counters.

## Comfortable and easy to use

In order to offer you the best ergonomics, the handles are shaped to distribute the weight you put on them evenly. The soft braking system ensures that minimal hand pressure is required, allowing users with reduced hand functions to easily lock the wheels in place.



## Accessories

	<b>Tray</b> Practical for resting items.  Part no, Fellow Classic: 08806 (image 1)  Part no, Fellow Supreme: 10667 (Image 2)		<b>Cane holder</b> Simplifies the transportation of a cane. Part no: 10603		<b>Oxygen holder</b> Part no, Fellow Classic: 08811.  Part no, Fellow Supreme: 10671
	<b>Basket with handle</b> Basket ideal for shopping.  Part no, Fellow Classic standard: 08805, Fellow Classic low: 88052. Fellow Classic narrow, wide and high: 88052. (Image 1)  Part no Fellow Supreme: 10666 (Image 2)		<b>Soft carry case</b> ideal for personal items. Part no: 08808. (Cannot be used on Fellow Classic low)		<b>Newspaper strap</b> Can only be used on Fellow Supreme. Part no. 10679.
	<b>Back rest</b> Gives support for the back when sitting. Part no: 10668		<b>Curb climber</b> Simplifies the transport over obstacles. Integrated in rollator for Fellow Supreme. Part no (Fellow Classic):10661		<b>One-hand brake</b> Allows control of brakes with one hand. Part no, Fellow Classic, right: 10677, left: 10678. Part no, Fellow Supreme, right: 10669, left: 10676.
	<b>Slow down brake</b> Prevents the rollator from moving too fast. Part no, right: 88121, Part no, left: 88122. Part no: 10625 (set)		<b>IV-stand</b> Can only be used with Fellow Classic. Part no:08814		

Technical Specifications									
Description	Part no.	Handle Height min	Handle Height max	Seat height	Width	Width between handles	Wheel size	Weight	Max user weight
Fellow Supreme	85400	73 cm / 28.7"	100 cm / 39.4"	590 mm / 23.2"	630 mm / 24.8"	440 mm / 17.3"	20 x 4 cm / 7.7"x1.7"	8.0 kg / 17.6 lbs	150 kg / 330.7 lbs
Fellow Classic std	88009	77 cm / 30.3"	95 cm / 37.4"	59 cm / 23.2"	63 cm / 24.8"	44 cm / 17.3"	16 x 4 cm / 7.7"x1.7"	8.0 kg / 17.6 lbs	150 kg / 330.7 lbs
Fellow Classic low	88109	68 cm / 26.8"	81 cm / 31.7"	51 cm / 19.9"	62 cm / 24.4"	44 cm / 17.3"		7.7 kg / 16.9 lbs	130 kg / 286.6 lbs
Fellow Classic wide	88309	77 cm / 30.3"	95 cm / 37.4"	60 cm / 23.6"	69 cm / 27.2"	50 cm / 19.7"		8.7 kg / 19.2 lbs	150 kg / 330.7 lbs
Fellow Classic narrow	88509	77 cm / 30.3"	95 cm / 37.4"	60 cm / 23.6"	60 cm / 23.4"	44 cm / 17.3"		8.1 kg / 17.8 lbs	130 kg / 286.6 lbs
Fellow Classic high	88609	82 cm / 32.3"	100 cm / 39.4"	64 cm / 25.2"	63 cm / 24.6"	44 cm / 17.3"		8.1 kg / 17.8 lbs	130 kg / 286.6 lbs



## RA-Rollators Trolley Walkers

Carl-Oskar RA	26
Rebel RA	26
Staffan RA	26
Carl-Oskar TW	28
Rebel TW	28
Staffan TW	28
Svea TW	28

# RA- Rollators



## Carl-Oskar RA

Carl-Oskar RA is designed for users requiring extra support or who find it difficult to maneuver an ordinary rollator. Ideal for those users recovering from stroke, the Carl Oskar RA features separate forearm supports and large ball-bearing mounted wheels that ensure walking is easy on uneven or rough surfaces.

The handles can be adjusted to many different heights and the supports are adjustable in angle, ensuring a comfortable fit. Folding is easy with the release handle that is integrated right in the seat.



## Rebel RA

A safe and reliable rollator with separate underarm supports, the Rebel RA is the ideal rollator for those pain-sensitive users suffering from arthritis, Parkinson's or many other neurological disorders.

Not only does it provide the extra support required to those with difficulty manoeuvring an ordinary rollator, its unique front fork design and wheels ensure easy and smooth walking. The Rebel RA is available with or without a seat and can be used indoors or outdoors.




## Staffan RA

Staffan RA is a highly maneuverable rollator for indoor use. It is only 59 cm wide, which makes it easy to move in tight spaces. Staffan RA also has an adjustable underarm support, providing a steady and ergonomic grip for the user.

Staffan RA has fixed rear wheel and brake for better stability and security. It has a maximum height of 128 cm (Staffan high) and can be folded for easier transport.



## Accessories

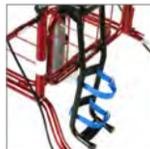	<b>Tray</b> Practical for resting light weight items. Part no: 10601. Not available for Staffan RA.		<b>One-hand brake</b> Allows the user to control brakes with one hand. Part no. Carl-Oskar RA, Rebel RA, Staffan RA low: 10653. Staffan RA, high: 10654		<b>Oxygen holder</b> Gives steady storage of the oxygen tube. Part no: 10615 (Staffan RA: 40605).
	<b>Slow down brake</b> Prevents the rollator from moving too fast. Part no. Carl-Oskar RA, Rebel RA: 10620. Not available for Staffan RA.		<b>Curb climber</b> Helps raise rollator onto a curb. Part no. Carl-Oskar RA, Rebel RA: 10612. Not available for Staffan RA.		<b>Basket with handle</b> Basket ideal for shopping. Part no: 10601 (Staffan RA: 40601).
	<b>Cane holder</b> Simplifies the transportation of a cane. Part no. Carl-Oskar RA, Rebel RA: 10603. Part no. Staffan RA: 40603		<b>Splash guard</b> Protects legs and clothes. Sold in pairs. Part no: 10608. Not available for Rebel RA and Staffan RA.		

## Technical Specifications

Description	Part no.	Handle Height Min	Handle Height Max.	Width	Wheel size	Weight	Max user weight
Carl-Oskar RA, Low 67cm	12315	95 cm / 37.4"	113 cm / 44.5"	67 cm / 26.4"	30 x 5 cm / 11.8"x1.8"	13.0 kg / 28.6 lbs	125 kg / 275.6 lbs
Carl-Oskar RA, Nor 67cm	12325	107 cm / 42.1"	126 cm / 49.6"	67 cm / 26.4"		13.2 kg / 29.1 lbs	
Carl-Oskar RA, Low 72 cm	12415	95 cm / 37.4"	113 cm / 44.5"	72 cm / 28.3"		13.1 kg / 28.9 lbs	
Carl-Oskar RA, Nor 72cm	12425	107 cm / 42.1"	126 cm / 49.6"	72 cm / 28.3"		13.3 kg / 29.3 lbs	
Rebel RA, Low 67cm	22311	86 cm / 33.8"	104 cm / 40.9"	67cm / 26.3"	20 x 4 cm / 7.9"x1.6"	11.3 kg / 24.9 lbs	125 kg / 275.6 lbs
Rebel RA, Nor 67cm	22321	104 cm / 40.9"	122 cm / 48.0"	67 cm / 26.4"		11.5 kg / 25.3 lbs	
Rebel RA, Low 72 cm	22411	86 cm / 33.8"	104 cm / 40.9"	72 cm / 28.3"		11.4 kg / 25.1 lbs	
Rebel RA, Nor 72 cm	22421	104 cm / 40.9"	122 cm / 48.0"	72 cm / 28.3"		11.6 kg / 25.6 lbs	
Staffan RA, Low 59cm	42221	86 cm / 33.8"	108 cm / 42.5"	59 cm / 23.2"	20 x 4 cm / 7.9"x1.6"	10.3 kg / 22.7 lbs	90 kg / 198.4 lbs
Staffan RA, High 59cm	42231	107 cm / 42.1"	128 cm / 50.4"	59 cm / 23.2"		10.6 kg / 23.4 lbs	

# Trolley Walkers

## Carl-Oskar Trolley Walker

With a wide, comfortable forearm support platform table, the Carl-Oskar Trolley Walker will easily assist those users requiring even more support than with an ordinary trolley walker. This comfortable and firm forearm table ensures those patients that are pain sensitive will be able to comfortably rest forearms while the front-facing brake handles will assist with maintaining better posture.



## Staffan Trolley Walker

Highly manoeuvrable, this trolley walker is designed for indoor use only with a narrow overall width of 62cm. Available in two heights, it folds easily for transport or storage.

The fixed rear wheels and brake add stability while the wide, firm platform table with raised side supports enhance comfort and safety.



## Rebel Trolley Walker

The Rebel Trolley Walker is a safe and reliable rollator with table pad that is designed for users requiring extra support or those who find it difficult to maneuver an ordinary walker. The unique front fork design and ball-bearing mounted wheels ensure very easy manoeuvrability – a primary feature for pain-sensitive users. For indoor or outdoor use.



## Svea Trolley Walker

This sturdy and safe trolley walker features a hydraulic gas spring that acts as a shock absorber while walking. Available in two widths and two heights, the Svea TW combines maximum weight bearing relief with the manoeuvrability of a rollator.

A stabilizer wheel for added stability is also available. The Svea is recommended for indoor use.



## Accessories



**Tray**  
Practical for resting light weight items. Part no. Carl-Oskar TW, Rebel TW: 10601. Not available for Staffan TW or Svea TW.



**One-hand brake**  
Allows the user to control brakes with one hand. Part no: Staffan TW, high: 10654. All other models: 10653.



**Oxygen holder**  
Gives steady storage of the oxygen tube, part no. 10615 (Staffan 40605). Not available for Svea TW.



**Slow down brake**  
Prevents the rollator from moving too fast. Part no. Carl-Oskar TW, Rebel TW: 10620. Not available for Staffan TW or Svea TW.



**Curb climber**  
Simplifies the transport over obstacles. Part no. 10612. Not available for Staffan TW or Svea TW.



**Basket with handle**  
Basket ideal for shopping. Part no: 10601 (Staffan TW 40603)



**Cane holder**  
Simplifies the transportation of a cane. Part no: 10603 (Staffan TW: 40603). Not available for Svea TW.



**Splash guard**  
Protects legs and clothes. Sold in pairs. Part no: 10608. Only for use with Carl-Oskar TW.

## Technical Specifications

Description	Part no.	Height Min	Height Max.	Width	Wheel size	Weight	Max user weight
Carl-Oskar TW, Low 67 cm	13315	95 cm / 37.4"	113 cm / 44.5"	67 cm / 26.4"	30 x 5 cm / 11.8"x1.8"	14.5 kg / 31.9 lbs	125 kg / 275.6 lbs
Carl-Oskar TW, Nor 67cm	13325	107 cm / 42.1"	126 cm / 49.6"	67 cm / 26.4"		14.7 kg / 32.4 lbs	
Carl-Oskar TW, Low 72 cm	13415	95 cm / 37.4"	113 cm / 44.5"	72 cm / 28.3"		14.6 kg / 32.2 lbs	
Carl-Oskar TW, Nor 72cm	13425	107 cm / 42.1"	126 cm / 49.6"	72 cm / 28.3"		14.8 kg / 32.6 lbs	
Rebel TW, Low 67cm	23311	86 cm / 33.8"	104 cm / 40.9"	67 cm / 26.4"	20 x 4 cm / 7.9"x1.6"	12.8 kg / 28.2 lbs	125 kg / 275.6 lbs
Rebel TW, Nor 67cm	23321	104 cm / 40.9"	122 cm / 48.0"	67 cm / 26.4"		13.0 kg / 28.7 lbs	
Rebel TW, Low 72 cm	23411	86 cm / 33.8"	104 m / 40.9"	72 cm / 28.3"		12.9 kg / 28.4 lbs	
Rebel TW, Nor 72 cm:	23421	104 cm / 40.9"	122 cm / 48.0"	72 cm / 28.3"		13.1 kg / 28.9 lbs	
Staffan TW, Low 62cm	43211	86 cm / 33.8"	108 cm / 42.5"	62 cm / 24.4"	20 x 4 cm / 7.9"x1.6"	11.7 kg / 25.8 lbs	90 kg / 198.4 lbs 125 kg / 275.6 lbs
Staffan TW, High 62cm	43221	107 cm / 42.1"	128 cm / 50.4"			11.9 kg / 26.2 lbs	
Svea TW Low 59 cm	53211	82 cm / 32.3"	107 cm / 42.1"	59 cm / 23.2"	20 x 4 cm / 7.9"x1.6"	14.9 kg / 32.8 lbs	150 kg / 330.7 lbs
Svea TW Nor 59 cm	53221	101 cm / 39.7"	126 cm / 49.6"	59 cm / 23.2"		15.1 kg / 33.3 lbs	
Svea TW Low 68 cm	53411	82 cm / 32.3"	107 cm / 42.1"	68 cm / 26.8"		15.1 kg / 33.3" lbs	
Svea TW Nor 68 cm	53421	101 cm / 39.7"	126 cm / 49.6"	68 cm / 26.8"		15.3 kg / 33.7" lbs	

## How do I choose the right rollator?

Consider the following before purchasing a rollator:

### 1. Which rollator should I use?

There are rollators that are good for every environment and others more suited for just indoor use. If using your rollator outside or on uneven surfaces, choose a rollator with large and stable wheels. If you mostly use it indoors, smaller wheels ensure walking is smooth. If you frequently transport it, think about weight and folding possibilities.

### 2. What size should it be?

The seat should feel roomy with sufficient distance between the handles so you do not have to squeeze into the seat. Make sure the width of the rollator fits through your doorways. A rollator that can partially fold while walking is ideal for narrow doorways. Your feet should rest comfortably on the ground when you sit down on your rollator.

### 3. How do I choose the right handle?

With weak or painful hands or wrists, anatomically shaped handles will spread the weight over a wider area of the palm thus reducing stress and making the use of the brake much easier.

### 4. How do I set the correct handle height?

When you stand straight, the handles should meet your wrist and arms should be slightly bent at the elbows [about 20 to 30 degree angle].

### 5. What wheels are appropriate?

Larger wheels facilitate walking on uneven ground while small wheels are better for indoor use.

### 6. What accessories do I need?

We offer a wide assortment of accessories to make your rollator individually yours! From deep baskets that are ideal for your groceries, to soft bags that can be left on the rollator while it is partially folded, to cane holders and trays – all our accessories will make your day to day outings easier and more enjoyable!



## How to measure for a rollator



### Measure the handle height

To measure for handle height, start by standing straight up with your shoes on. Bend your arm to approximately a 20° to 30° angle. Have someone measure from the crease of your wrist to the ground. This measurement should be the approximate height of the handles.



### Measure the seat height

Sit in a chair that you are comfortable in and from which you can easily rise from. Measure the floor to seat height and this will give you an approximate seat height in a rollator.



### Measure your doorways

Find the smallest place you will need to fit through. Most often, that is a doorway or hallway. Ensure that your rollator's overall width will go through this space. A rollator that can partially fold while walking may also be used for those smaller spaces.



Distributor:



**HUMAN CARE**  
 Your Life. Your Way.

**AUSTRALIA**

Unit 4, 49 Bertie Street  
 Port Melbourne, VIC 3207  
 Phone : +61 1300 132 736  
 Fax : + 61 (03) 9645 4954

info.au@humancaregroup.com

**CANADA**

10-155 Colonnade Road  
 Ottawa, ON K2E 7K1  
 Phone: 613.723.6734  
 Fax: 613.723.1058

info.ca@humancaregroup.com

**SWEDEN**

Årstaängsvägen 21B  
 117 43 Stockholm  
 Phone: +46 8 665 35 00  
 Fax: +46 8 665 35 10

info.se@humancaregroup.com

**NETHERLANDS**

Elspeterweg 124  
 8076 PA Vierhouten  
 Phone: 0577-412171  
 Fax: 0577-412170

info.nl@humancaregroup.com

**UNITED STATES**

4210 S. Industrial Drive, Suite 160  
 Austin, TX 78744  
 Phone: 512.476.7199  
 Fax: 512.476.7190

info.us@humancaregroup.com